

PRODUCT SPECIFICATION

Product PN 10496114

Mod. 984

Description M-FC with Rosolic Acid Broth - Ampoules

Rev. 02

M-FC with Rosolic Acid Broth - Ampoules



PRODUCT DESCRIPTION	AMP M-FC w/ Rosolic Acid 2ml, (50pk) M-FC with Rosolic Acid acts and functions in the same way as M-FC Broth (which allows the development of fecal coliforms at elevated 44.5°C temperature). Rosolic acid inhibits bacterial growth in general, except for fecal coliforms. Bile salts inhibit non-enteric bacteria. Aniline blue indicates the ability of fecal coliforms to ferment lactose to acid by-products that causes a pH drop in the medium. Other organisms form grey to cream colonies.
MANUFACTURER NAME	GVS North America 63 Community Drive Sanford, Me 04073 Phone: +1.866.736.1250 eMail: CustomerCareGVSLs@gvs.com - Website: www.gvs.com
INTENDED USE / APPLICATION	M-FC with rosolic acid allows the development of fecal coliforms at elevated temperatures (44.5°C).
MATERIALS	Media: m-FC with Rosolic Acid Broth Color: Reddish-purple nutrient liquid Ampoules: polypropylene Regulatory Documentation Required: ☐ Biocompatibility according ISO 10993-1 ☐ IMDS ☒ RoHS, directive 2002/95/CE ☐ BSE/TSE, directive 2003/32/ ☐ DEHP plasticizer Free and latex free ☐ Aging ☒ Reach 1907/2006/CE (hazardous substances regulation) ☐ Dir. 67/548/CE and Reg. 1272/2008/CE (medical sector dangerous substances) ☒ Conflict minerals

PRODUCT SPECIFICATION

Product PN **10496114**

Mod. 984

Description **M-FC with Rosolic Acid Broth - Ampoules**

Rev. 02

PRODUCT CHARACTERISTIC	Reddish-purple liquid media dispensed in plastic ampoules. Physical: Ampoule Dimensions: 57.9 mm x 8.4mm (at cap)/22 mm (at seal) Volume/ampoule: 2 - 2.4 ml Incubation Temperature Range: 44.5°C ± 2°C Storage temperature Range: From 2°C to 8°C Biological: Product is sold sterilized. Toxicity: See SDS for further information. Media Formulation: Approximate formula* per liter of water <table border="1"> <tbody> <tr> <td>Tryptose</td><td>10.0 gm</td></tr> <tr> <td>Peptone</td><td>5.0 gm</td></tr> <tr> <td>Yeast Extract</td><td>3.0 gm</td></tr> <tr> <td>Sodium Chloride</td><td>5.0 gm</td></tr> <tr> <td>Lactose</td><td>12.5 gm</td></tr> <tr> <td>Bile Salts</td><td>1.5 gm</td></tr> <tr> <td>Aniline Blue</td><td>0.1 gm</td></tr> <tr> <td>Rosolic Acid, 1%</td><td>10 ml</td></tr> </tbody> </table> <i>* Adjusted and/or supplemented as required to meet performance requirements pH 7.4 ± 0.2</i>	Tryptose	10.0 gm	Peptone	5.0 gm	Yeast Extract	3.0 gm	Sodium Chloride	5.0 gm	Lactose	12.5 gm	Bile Salts	1.5 gm	Aniline Blue	0.1 gm	Rosolic Acid, 1%	10 ml
Tryptose	10.0 gm																
Peptone	5.0 gm																
Yeast Extract	3.0 gm																
Sodium Chloride	5.0 gm																
Lactose	12.5 gm																
Bile Salts	1.5 gm																
Aniline Blue	0.1 gm																
Rosolic Acid, 1%	10 ml																
INSTRUCTIONS WARNINGS	Product is temperature sensitive. Store at recommended temperature.																
PRODUCT SHELF LIFE	9 months when stored properly																
STERILIZATION	Product is sold sterilized.																

PRODUCT SPECIFICATION

Product PN **10496114**

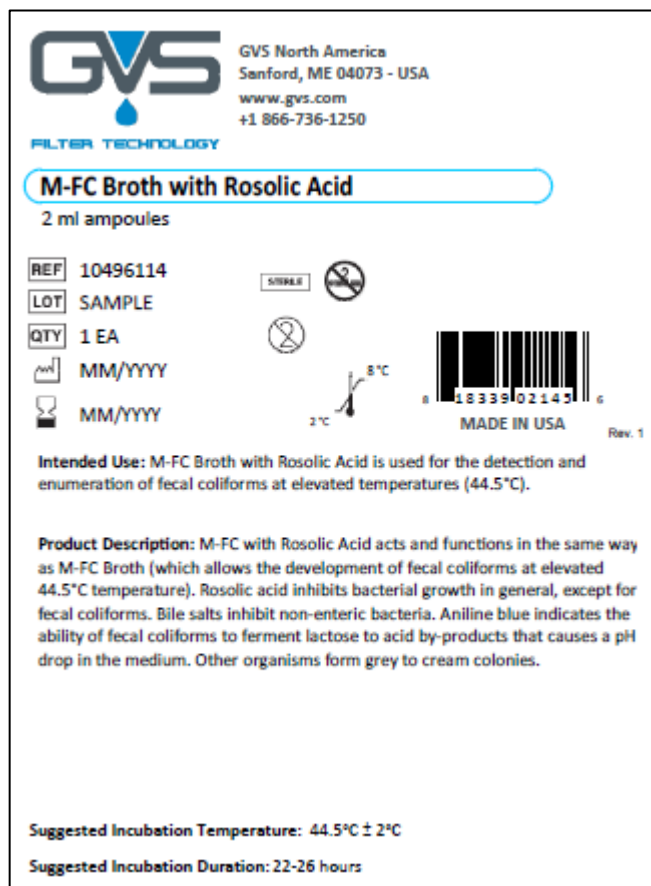
Mod. 984

Description **M-FC with Rosolic Acid Broth - Ampoules**

Rev. 02

PACKAGING AND LABELLING

Primary: 50 ampoules are placed in a labeled foil bag and sealed using a heat sealer.



CERTIFICATE OF COMPLIANCE

With each shipment, GVS NA will provide the CoC with item number, lot number and expiration date. The Quality management system is in compliance with ISO 9001.

VISUAL REQUIREMENTS

Visual acceptance requirements apply when inspected under below conditions:

Magnification: none

Light type: normal room lighting

Acceptance Requirement			Sampling Plan
1	Color – reddish purple	all	100%
2	Ampoule appearance – no misshapen or leaking ampoules	none	

PRODUCT SPECIFICATION

Product PN 10496114

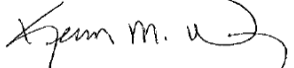
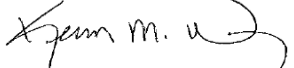
Mod. 984

Description M-FC with Rosolic Acid Broth - Ampoules

Rev. 02

PERFORMANCE REQUIREMENTS	Acceptance Requirement		Sampling Plan
	Parameter	Specification	ISO 11133:2014
	1 Recovery (Measured against control)	≥ 85%	ISO 11133:2014
	2 Positive Test Organism	<i>E.coli</i> ATCC 25922	ISO 11133:2014
	3 Negative Test Organism	<i>Ps. aeruginosa</i> ATCC 10145	ISO 11133:2014
	4 pH at 25°C	7.4 ± 0.2	ISO 11133:2014
	5 Sterility	Zero Growth	ISO 11133:2014
	6 Media Color	Reddish-purple, Clear	ISO 11133:2014
	7 Colony Color	Fecal coliforms are blue in color, all others are gray to cream.	ISO 11133:2014
	8 Incubation	44.5°C for 22-26 Hours	ISO 11133:2014
This material specification describes the properties of product above indicated. This document contains general requirements, material description, drawing references, defect specification, biological material requirements.			

REVISIONS AND APPROVALS:

DATE	REV.	REASON FOR CHANGE	ISSUED AND CONTROLLED BY: (name /function and signature)	APPROVED BY: (name /function and signature)
11/7/18	0	Initial release	Debra English/Manager Biochemistry Applications 	Kevin Wrigley/Director, Quality 
7/28/2021	1	Update description, label, ampoule dimension & add media formulation	Debra English/Manager Biochemistry Applications 	Kevin Wrigley/Director, Quality 

Customer Approval:

PRODUCT SPECIFICATION

Product PN **10496114**

Mod. 984

Description **M-FC with Rosolic Acid Broth - Ampoules**

Rev. 02

We accept this material specification as a part of the agreed terms of delivery

Company name _____

Approved by: _____
(Name, Function) (Signature)

Date _____
(Company stamp)

Please send back this document signed for approval. If we will not receive this specification signed, we consider the first order placed as implicit approval.